


Piko[®] Plate Illuminator

Simplify
sample
loading

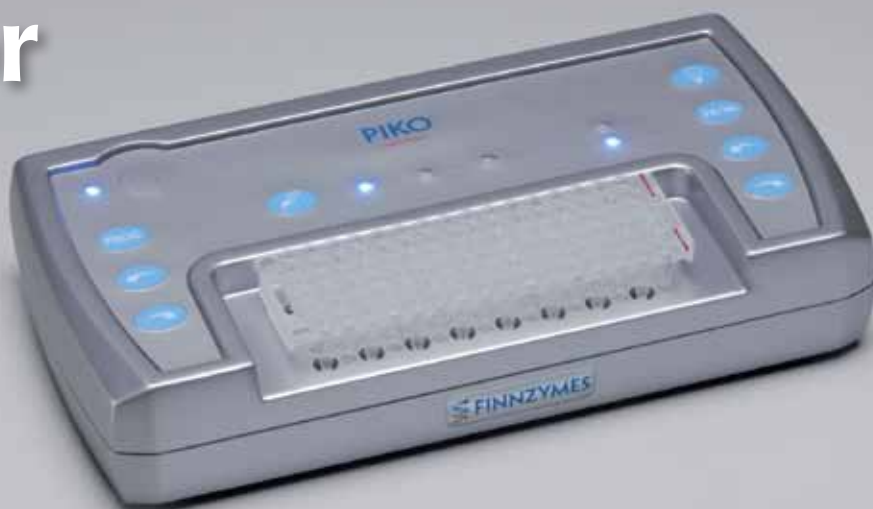

Pipetting aid for Piko[®] PCR Plates

Finnzymes' Piko[®] Plate Illuminator is a state-of-the-art plate rack that provides a simple way to track the loading of reactants by illuminating the target well(s) from below with white LED lights. With 10 different pre-programmed loading patterns, this device can fulfill just about any protocol requirements. It accepts both 24-well and 96-well Piko[®] PCR Plates and can be used with all standard single-channel, 8-channel or 16-channel pipettors.

Using Piko[®] Plate Illuminator is easy:

1. Select loading program
2. Place your plate in the rack
3. Start dispensing



Features and benefits

- **Compatibility** – works with all standard single channel, 8-channel and 16-channel pipettors
- **Ease of use** – compact and simple instrument with pre-programmed loading patterns
- **2 in 1 design** – both 24-well and 96-well Piko PCR Plates fit into the rack

Specifications	
Size	Width 16 cm, depth 9 cm, height 3 cm
Weight	0.1 kg
Vessel types	24-well and 96-well Piko PCR Plates
Programs	10 preset and user-selectable loading patterns
Power requirements	100–240 V input power

Ordering information

Product code	Description	Quantity
PIP2496	Piko [®] Plate Illuminator	1

www.finnzymes.com

Piko[®] is a trademark or registered trademark of Finnzymes Oy or its affiliates.